X Pro Elite

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Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Prague 7, www.alza.cz

Important Safety Notice


This product is recommended for use by persons over 12 years of age.

The product can only be used with Oclean replacement heads. We will not be responsible for any consequences resulting from the use of non-certified third-party brush heads.

This product is for personal cleaning only and is not recommended for medical purpose or multi-user use.

This product is for indoor use only. Do not place the toothbrush under direct sunlight for a prolonged period. It is recommended to be placed in the dry well-ventilated areas with normal temperature.

Please keep the body dry before charging, dry it with a soft cloth, and clean the charging base with a dry cloth regularly to prevent the stain from building up.

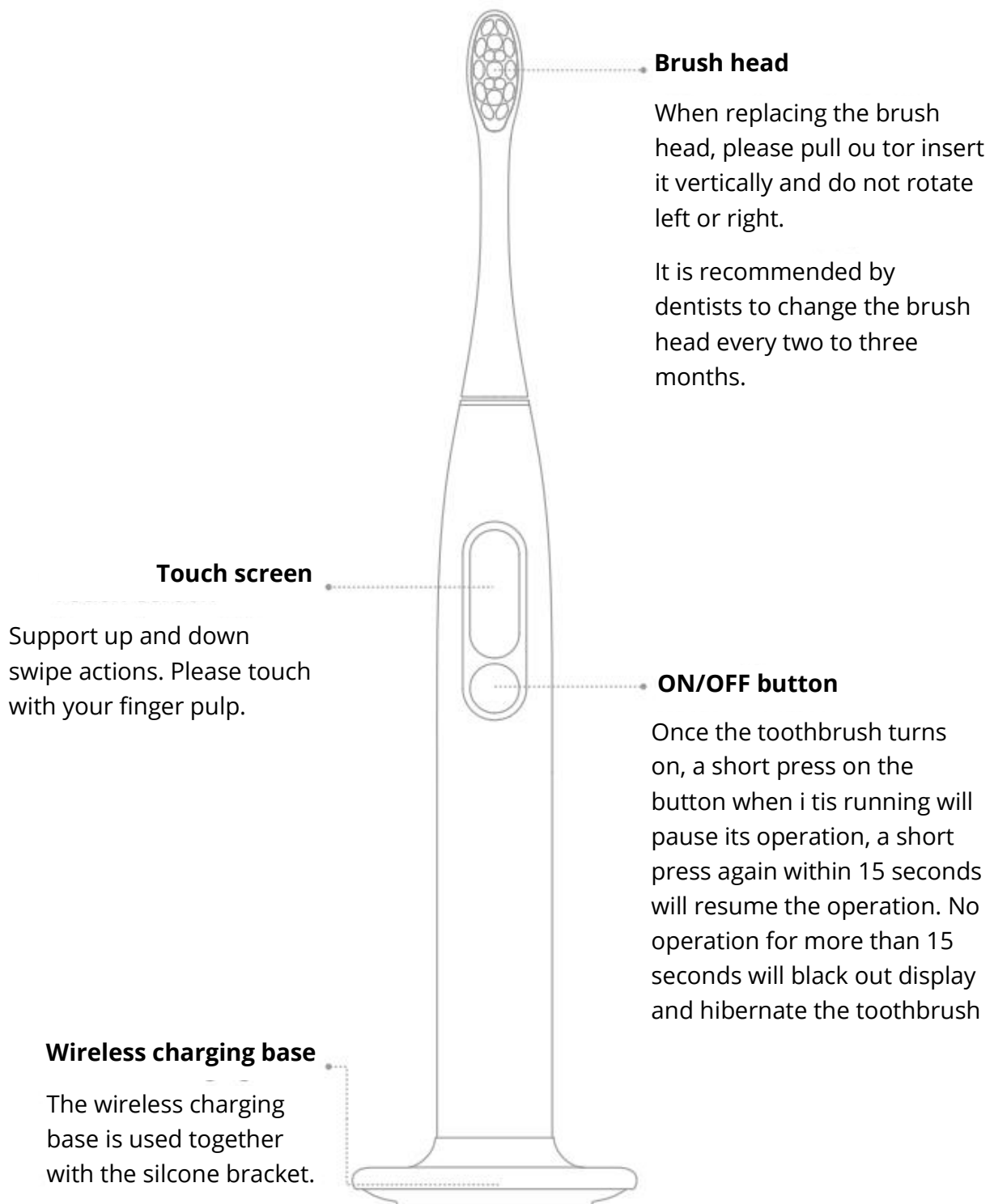Use for the first time

Tap the on/off button to turn on the toothbrush and enable Bluetooth. The toothbrush will then remain in standby mode. Bluetooth will also be in always-on mode.

Brushing program

Without any settings, the toothbrush will run for 2 minutes by default following “Standard Daily Cleaning Program”, where the toothbrush will alert the user of the need to change the brushing area by changing the vibration pattern every 30 seconds and will automatically stop running after 2 minutes of operation.

Charging the toothbrush

Dry the body before charging and then place it on the wireless charging base for charging. You can determine the charging status with the indicator light on the wireless charging base and its battery level with the battery percentage shown on the Oclean electric toothbrush display. Please do not leave the toothbrush on the wireless charging base for a prolonged period for repeated charging. After fully charged, please remove the toothbrush from the wireless charging base in time and recharge it when the toothbrush alerts your its low battery status.

Install the app

After scanning the QR code to download the Oclean APP, please search and connect your toothbrush in the app. Note: Attempt to connect the toothbrush in “System Settings – Bluetooth Settings” in your phone will not work.

Wireless charging base

Lay flat when charging and do not place mental other than toothbrush on the charging base. During charging, please mount the wireless charging protective sleeve on the charging base and then place the toothbrush in the slot of the wireless charging base.

Wall mount

The toothbrush comes with a wall mount for space-saving storage.

Note

Installation of the bottom of the wall mount

When placing the toothbrush on the wall mount, make sure it is facing up.

If you need more detailed product usage information, please follow the Oclean Official Account – www.facebook.com/ocleanofficial.

This is a smart sonic electric toothbrush product. Please read the user manual carefully to get started quickly. It can be used with the matching app. See below to learn more about brushing teeth with this toothbrush.

Attention should be paid for those with the following conditions or oral sensitivities

- Those who have recently undergone periodontal surgery or are undergoing treatments.
- Individuals with other symptoms in the mouth including those with unstable dentures or have persistent heavy bleeding from the gums.
- Other individuals with sensitive and fragile oral issues. This product is not recommended for use with above group or consult a profession dentist before use.
- This product is not suitable for people with disabilities, mental illness, or lack of relevant experience and knowledge [including children (0-12 years old)]. Unless used under the supervision of their guardians.

Instruction for use

Turning the product on for the first time

To use the toothbrush for the first time, please short press “On” button, and it will turn on the toothbrush, and Bluetooth will also be in always-on mode. You can download and install the “Oclean App” and connect and set the toothbrush in the app. App installation and toothbrush setup.

Scan the QR code to complete the app installation, which will allow you to set various functions for your toothbrush, such as: Turning on/off the “charging brushing alert after 30 seconds”, turning on/off the “Automatic Frequency Reduction Tooth Guard” function and etc. In addition, you can also customize your own brushing program and push the customized program directly to your toothbrush.

Touch the screen to setup the toothbrush

The toothbrush setup can also be made via the touch screen without connection to the app. You can choose between four modes: Cleaning, Message, Whitening and Sensitive as well as select the duration and intensity. After receiving the pushed brushing program from the app, if any of the mode, intensity or duration changes are made via the touch screen, then toothbrush will no longer perform the previous app-pushed brushing programs.

Timer and brushing area cleaning alert

The default duration of the brushing program is 2 minutes, with a momentary changing in vibration status every 30 seconds to alert the user of the need to change the brushing area.

Tooth brushing report

After the brushing is done, the product display will show your brushing time, brushing scores, and tips for blind spot of tooth brushing for the brushing you have just completed.

Brushing Mode Setup

Go to settings by swiping up and down the display or swiping left and right of the screen edge – e.g.: The “Brushing Mode” setting. The toothbrush will automatically remember the settings once the setup is completed, and the subsequent operation will follow the new settings.

For the intensity settings, swipe the screen up and down to change the intensity level to “plus 1” or “minus 1”, and swipe the screen edge up and down to change the intensity level to “plus 5” or “minus 5”

Swipe up and down

Press ON/OFF button confirm the settings

Adjusting the time

After connecting with the app and synchronizing the information of the paired toothbrush, the time on the toothbrush will be automatically updated to local time.

- Do NOT disassemble, strike, extrude or throw the device into a fire.
- Do NOT use the toothbrush if its battery is bulging or as expanded.
- Do NOT subject the battery to a high temperature environment.
- Do NOT use the battery if it has been subjected to moisture.
- Do NOT place metal other than toothbrush on the charging base.
- Before appliance is discarded, its battery must be removed from it.
- Before the battery is removed, the appliance must be powered off.
- The battery shall be disposed of safely.

- New Message Icon
The new message icon will appear on the display when there is new activity message, software update notification or another alert information.
- Check the Brushing Records in the App
The toothbrush must be paired with and connect to the app in order to receive the brushing data from the brushing. As the toothbrushing comes with limited memory size for data storage, the brushing data can usually be stored for no more than 7 days.

- Low Battery Alert
When the battery percentage is below 20%, an alert will pop on the display of your toothbrush from time to time to indicate that the battery is running low, and you need to charge it soon. When the battery is seriously low, and the percentage is lower than 1%, your toothbrush will be off, and the toothbrush's motor will vibrate to alert you to charge the battery, and the screen will not be turned on either.


How to Charge

Please place the designated wireless charging base of the product on a smooth surface and connect the charging base to the plug. When power is applied, the indicator light on the base will flash. After mounting the silicone protective sleeve on the surface of the wireless charging base, place the toothbrush handle in the slot of the wireless charging base. If not charged, the indicator light will not lit. When there is a foreign body on the wireless charger or when it is not being charged properly, the indicator light will flash quickly. When the toothbrush is being charged properly, the indicator light will turn on a breathing-like display effect and the battery will generally be fully charged in around 3.5 hours.

Important Charging Notice

This product adopts wireless charging mechanism, it is recommended to place this product on the wireless charger only when it needs to be charged, when it does not need to be charged, please remove the toothbrush from the charging base. Leaving the toothbrush on the powered wireless charging base for a prolonged period of time may lead to premature depletion of the build-in battery affecting the usage experience.

Recommended tooth brushing method

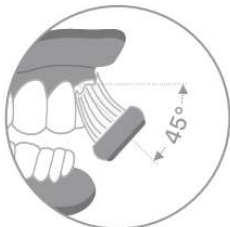Move the toothbrush slowly at 45° between the brush and gingiva line. There is no need to manually brush up and down during this process.

Hold the toothbrush handle gently, move slowly in a circular motion, and please do not apply excessive pressure.

Please gargle with clear water after brushing.

Cleaning and Maintenance

Rinse with tap water and air dry, and do not pull the bristles of the brush head.

Product specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Model	Y2087
Version	Pro Elize
Emission power	0 dBm
IP Rating	IPX7
Power of the toothbrush body	1.5W
Charging	Oclean USB charging base
Input	5V $\equiv$ 0.3A
Operating time*	20 days (level 32 of intensity) – 40 days (level 1 of intensity)
Battery Capacity	800mAh
Charging time	3.5h
Implementation standard	GB 4706.1-2005 GB 4706.59-2008
Bluetooth	Bluetooth 5.0 above

* Under the cleaning mode, it is calculated at 2 minutes/time, 2 times/day working condition.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

CIN: 27082440

Subject of the declaration:

Title: Smart Sonic Electric Toothbrush

Model / Type: Y2087

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/35/EU


Directive No. 2014/53/EU

Directive No. 2011/65/EU as amended 2015/863/EU

WEEE

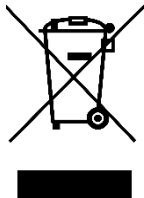This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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